

# ***Manuals for the EVAN Toolbox No. 1***

V2.3 - November 2017

**Paul O'Higgins<sup>1</sup>, Gerhard W. Weber<sup>2</sup>,  
Cinzia Fornai<sup>2</sup>, Fred L. Bookstein<sup>2</sup>**

<sup>1</sup> University of York

<sup>2</sup> University of Vienna

## **ET Templand (ET 1.72)**

**Creating landmarks and sliding semilandmarks on curves and surfaces**

<b>1. The EVAN Toolbox (ET)</b>	<b>2</b>
<b>2. Landmarks and semilandmarks</b>	<b>2</b>
<b>3. The example data</b>	<b>2</b>
<b>4. File Types</b>	<b>5</b>
<b>5. The principle workflow</b>	<b>6</b>
<b>6. Starting Templand</b>	<b>6</b>
<b>7. ET Templand Graphical User Interface (GUI)</b>	<b>7</b>
<b>8. Loading a Surface or a Curve</b>	<b>8</b>
<b>9. Placing a landmark</b>	<b>10</b>
<b>10. Adding semilandmarks</b>	<b>10</b>
<b>11. Save and export a form file</b>	<b>11</b>
<b>12. Creating a tableau file</b>	<b>11</b>
<b>13. Working with the example data</b>	<b>12</b>

The EVAN Toolbox is software designed by Roger Phillips of the University of Hull, Paul O'Higgins of the University of York, and William D. K. Green and Fred Bookstein of the University of Vienna. Its principal systems programmers were Helgi Gunnarsson, Ramy Gowigati, Youssef Shady George-Nashed, Vincent Dalge, and Oualid Ben Ali.



## **1. The EVAN Toolbox (ET)**

The EVAN Toolbox (ET) is a software package developed by the European Virtual Anthropology Network - EVAN ([www.evan.at](http://www.evan.at)) and the EVAN-Society **to facilitate form and shape analysis of objects featuring a complex geometry**. It uses Geometric Morphometrics (GM) which includes methods such as General Procrustes Analysis, Principal Component Analysis, Thin-Plate Spline Warping or Partial Least Squares Analysis. These features are all embedded in the **ET Core**.

In addition, the EVAN Toolbox supports **data acquisition**, i.e. to locate landmarks and semilandmarks on surfaces and curves. This part is realized in **ET Templand**, a node that can be started from the ET user interface. This manual describes how to use this Templand node.

## **2. Landmarks and semilandmarks**

**ET Templand** is a tool that allows landmarking and semilandmarking of specimens. While landmarks (LM) characterize biologically homologous points on a specimen (such as the crossing of the sagittal and coronal sutures or the most extreme point on the mastoid process), sliding **semilandmarks** (sLM) are placed along **curves** (e.g. on ridges like the orbital rim) or on **surfaces** (e.g., on the parietal bone where not much real landmarks are) (see the six types of landmarks in Weber & Bookstein, Springer Verlag, 2011). Sliding semilandmarks are considered to be geometrically homologous (Gunz & Mitteroecker 2013; Gunz et al. 2005). The placing of semilandmarks requires a sample. Semilandmarks are based on the Thin Plate Spline (TPS) function (Bookstein 1991) and allowed to slide in relation to the landmark configurations in the sample. This sliding involves different degrees of freedom:


- a) Semilandmarks on curves have 1 degree of freedom to slide  
(along the defined curve)
- b) Semilandmarks on surfaces have 2 degrees of freedom to slide  
(on the defined surface)
- c) Points that are declared "missing" have 3 degrees of freedom to slide  
(within the whole data volume)

The ET Templand node offers procedures to identify classic landmarks and also semilandmarks mentioned under point a) and b) above. To start the process, a **template** is defined. At start this will be any specimen, preferably a specimen **near the mean shape**. All LM and sLM are identified on this specimen. Then, we have to identify LM on a second specimen, the **first target** specimen. The sLM from the template will be projected to the target and then slid. This process is iterative until no differences are appreciable (usually, a few rounds of project-slide will do). The next specimen, the **second target**, will undergo the same procedure, and so on. Once all specimens in the sample have been done, the coordinates of LM and sLM are used for a "Procrustes Superimposition" (which you know from the standard features of ET Core) and the **template can be replaced by the resulting average configuration** (currently by hand and not yet automatically, but will be implemented for the next version). All specimens are then slid against this average again. These final results can be used for **further analyses** of form and shape. For that, you will have to export them back to the **ET Core** where all the nodes for Procrustes Superimposition, Principal Component Analysis or Partial Least Squares analysis are available.

## **3. The example data**

We will use an example from dental anthropology to demonstrate the principle workflow and the individual working steps in this manual. Our example data involves of the dentin body of **human upper molars**. The

surfaces have been extracted from  $\mu$ CT images by segmenting the enamel and the dentin (which can be done e.g. in Amira, Avizo, TurtleSeg, ImageJ or similar software packages), creating triangulated surfaces for the

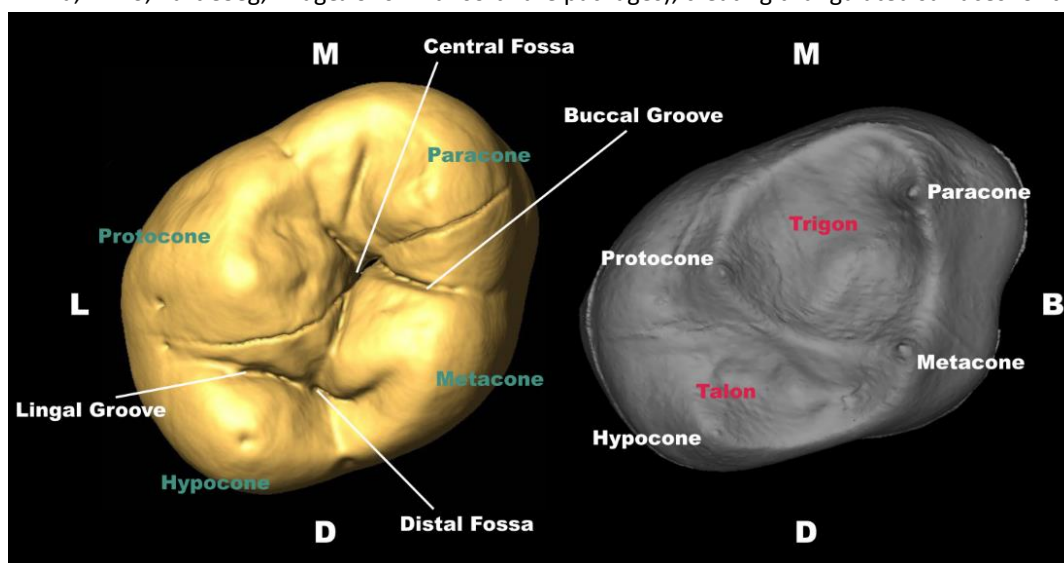

**dentin**, and saving these surfaces in **.obj format** (ET Templand will only use OBJ format for surfaces, ET Core can also handle other formats such as STL). For your orientation, we provide an annotated picture here showing the enamel surface (left - yellow) and the dentin surface (right -grey) of a human upper left second molar ( $M^2$ ) from occlusal view. M – mesial, D – distal, B – buccal, L – lingual.

The cusps and ridges on the dentin surface are easily identifiable. There are **10 classic landmarks** (yellow spheres in the figure below) which we will use in our data set. Each of them comes with a specific name and a definition where they are located. Choosing these landmarks on the surface of the dentin crown, we took care that they were **well distributed across the whole structure** because they establish a framework for sliding semilandmarks later.

Classic Landmarks		
No.	Landmark	Location
1	Cervix_Bucc	the midpoint of the cervical margin buccally
2	Cervix_Dist	the midpoint of the cervical margin distally
3	Cervix_Ling	the midpoint of the cervical margin lingually
4	Cervix_Mesi	the midpoint of the cervical margin mesially
5	Hypocone	the tip of the hypocone
6	Metacone	the tip of the metacone
7	Paracone	the tip of the paracone
8	Protocone	the tip of the protocone
9	Talon_Bott	the deepest point in the talon
10	Trigon_Bott	the deepest point in the trigon

The dentin crown also features some nicely expressed **ridges** which are ideal for curve semilandmarks (comparable to muscle markings on long bones or the temporal line on the cranium, if working with bones).

One of these ridges runs between the Protocone, the Paracone, and the Metacone. We call it **“Curve\_Trigon”** in our example and it consists of **15 curve semilandmarks** (blue spheres in the figure below).


The other curve can be traced by **5 curve semilandmarks** from the Metacone to the Hypocone and is called “**Curve\_Talon**” (green spheres).

Finally, we include **8 surface semilandmarks** (red spheres) within the trigon around the deepest point to capture the geometry of this depression “**Surface\_Trigon**”.

**Curve semilandmarks** have to be **affiliated to a predefined curve** which needs to be digitised in advance by following an actual ridges on the surface. **Surface semilandmarks** have to be **affiliated to a surface** which is the entire surface of the dentin crown in this case.

Curve_Trigon Semilandmarks		
No.	Semilandmark	Location
1	SM00001	start near the Metacone
2	SM00002	continue in the direction of the Protocone
3	SM00003	continue in the direction of the Protocone
4	SM00004	continue in the direction of the Protocone
5	SM00005	almost reached the Protocone
6	SM00006	start near the Protocone (on its other side)
7	SM00007	continue in the direction of the Paracone
8	SM00008	continue in the direction of the Paracone
9	SM00009	continue in the direction of the Paracone
10	SM00010	almost reached the Paracone
11	SM00011	start near the Paracone (on its other side)
12	SM00012	continue in the direction of the Metacone
13	SM00013	continue in the direction of the Metacone
14	SM00014	continue in the direction of the Metacone
15	SM00015	almost reached the Metacone
Curve_Talon Semilandmarks		
No.	Semilandmark	Location
1	SM00001	start near the Metacone
2	SM00002	continue in the direction of the Hypocone
3	SM00003	continue in the direction of the Hypocone
4	SM00004	continue in the direction of the Hypocone
5	SM00005	almost reached the Hypocone
Surface_Trigon Semilandmarks		
No.	Semilandmark	Location
1	SM00001	start near the Metacone
2	SM00002	continue counterclockwise
3	SM00003	continue counterclockwise
4	SM00004	continue counterclockwise
5	SM00005	continue counterclockwise
6	SM00006	continue counterclockwise
7	SM00007	continue counterclockwise
8	SM00008	continue counterclockwise


Please note: **For any study using landmarks, it is very crucial to make sure that the order of digitized landmarks is the same for each specimen.**

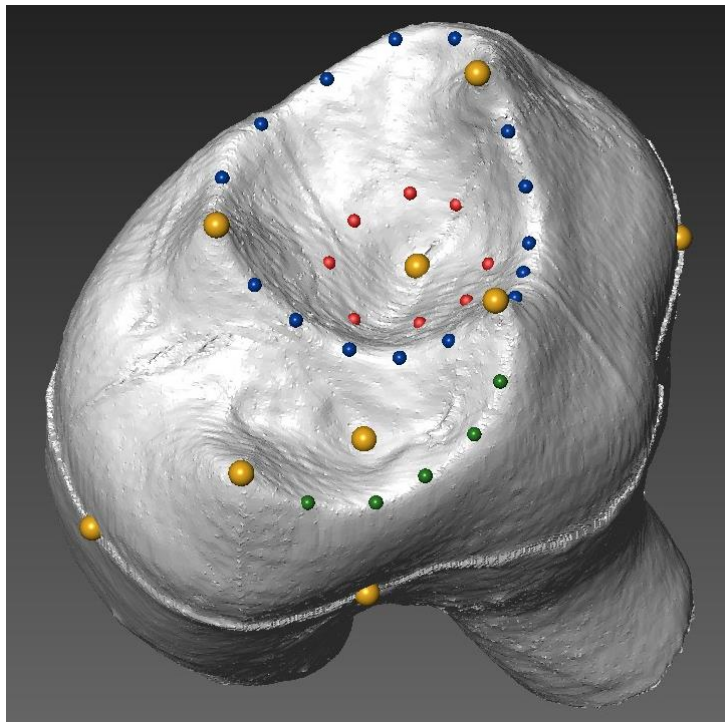

#### 4. File Types

Templand uses two types of file, one is a **form file** and the other is a **tableau file** (see figure below).

- **Surface files (.obj)** contain the 3D surface of objects which have been created in other software in .obj format (v-vertices, e-edges, f-faces). They can be edited in Wordpad if the format is distorted.
- **Form files (.frm)** contain **one individual specimen**. This will include the surface files (.obj), predefined curves (.obj), identified LMs and sLMs on surfaces and curves, and other information (e.g. about symmetry) related to the landmarking operation. Each form file can be thought of as storing either the full set of landmarks or an intermediate stage in landmarking.
- **Tableau files (.ltb)** contain information **on the sample**. They connect sequences of form files allowing, for example, the export of all landmark configurations at once or an animation or video of the landmarking sequence to be produced and displayed. Form files can be “appended”.

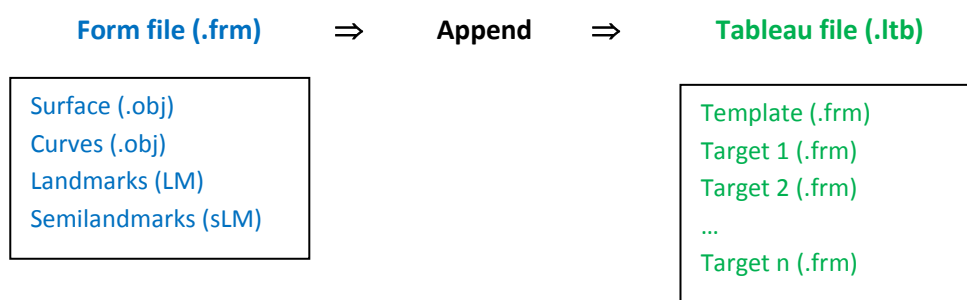

## **5. The principle workflow**

ET Templand provides the freedom of digitizing classic landmarks and preparing curves in any software and then import them into a form-file. However, Templand offers these functions as well. In our suggested workflow here we assume that all steps are carried out in Templand.


Once you have the surface files of all objects in the sample ready in .obj format, you can start digitizing classic landmarks. If you will have curve semilandmarks as well, it is however better to start building those curves first. For our example data, the workflow looks as follows:

- **Preparation of the sample**
  - Create surface files for each specimen as .obj in an appropriate software package
  - Create curves for each specimen and save them as .obj
  - Locate classic landmarks on each specimen, import curves, and save each form as .frm
  - Decide which specimen will be used as template (in the example data this is tooth "152")
  - Locate semilandmarks on the template specimen and save the .frm
- **Creating sliding semilandmarks**
  - Load the template form
  - Save a tableau file as .ltb
  - Append the template form as the first specimen
  - Change view to 4-window view ("Template")
  - Load the template form in the upper window
  - Reset the lower window and load the first target form (tooth "697")
  - Warp all unmatched landmarks from the template to the target
  - Append the first target to the tableau file
  - Reset the lower window and load the second target form (tooth "S4")
  - Warp all unmatched landmarks from the template to the target
  - Append the second target to the tableau file
  - ... continue for the whole sample ...
  - Save the tableau file
- **Analyzing the data**
  - Export landmarks for all specimens or import tableau file in ET Core
  - Perform analyses of your sample in ET Core


## **6. Starting Templand**

ET Templand is one of the available nodes in ET Core. To access Templand:

1. Start the EVAN toolbox
2. Go to the file menu
3. Select "New" to begin a new untitled network.
4. Click on the Templand node (tab "Data" of the left hand menu) just once with left mouse button
5. Then click on the work space just once with left mouse button
6. Double click on the Templand node (that appears now centre screen)
7. In the new window (Templand GUI) that appears click again "File" > "New"

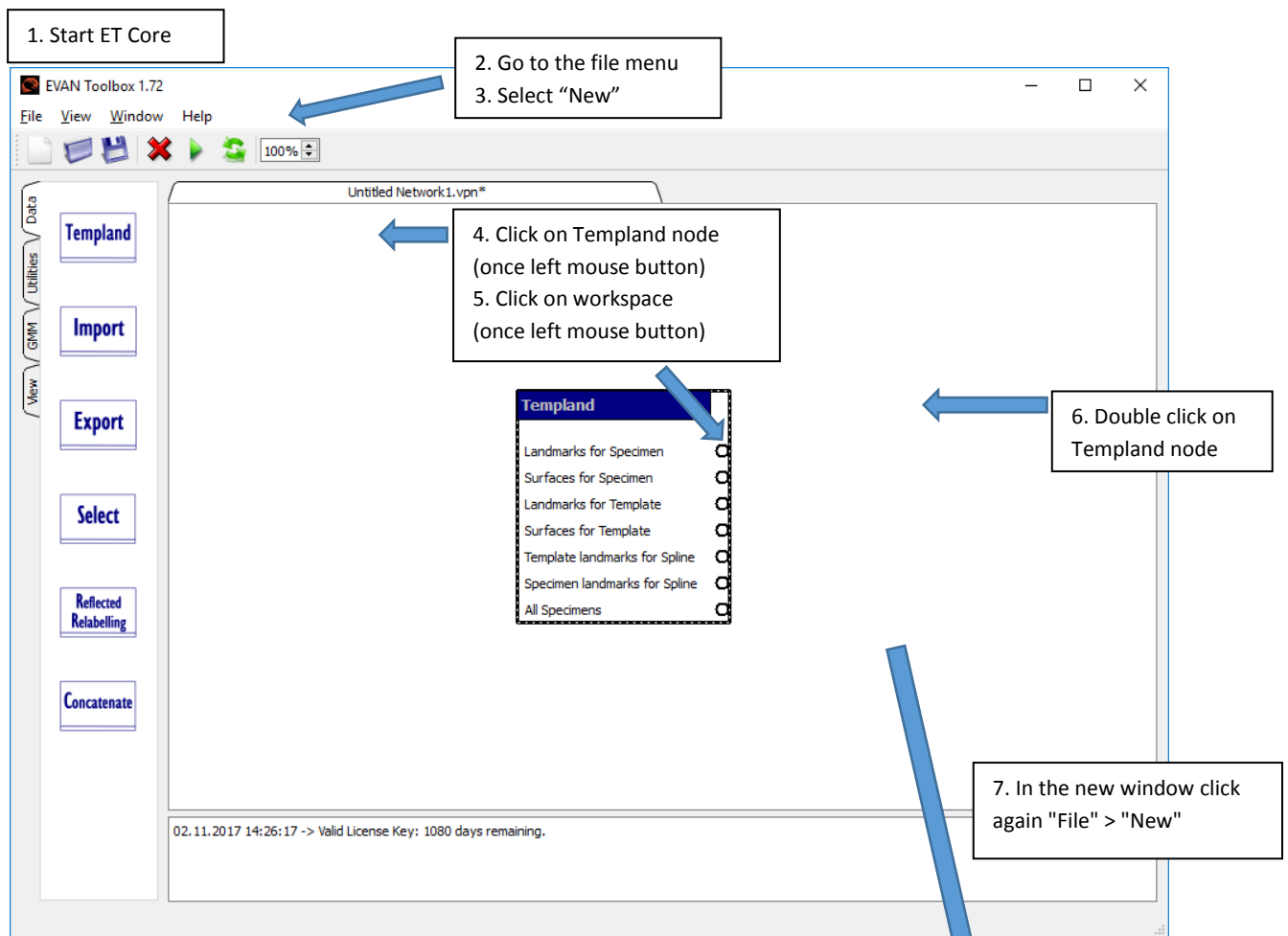

## 7. ET Templand Graphical User Interface (GUI)

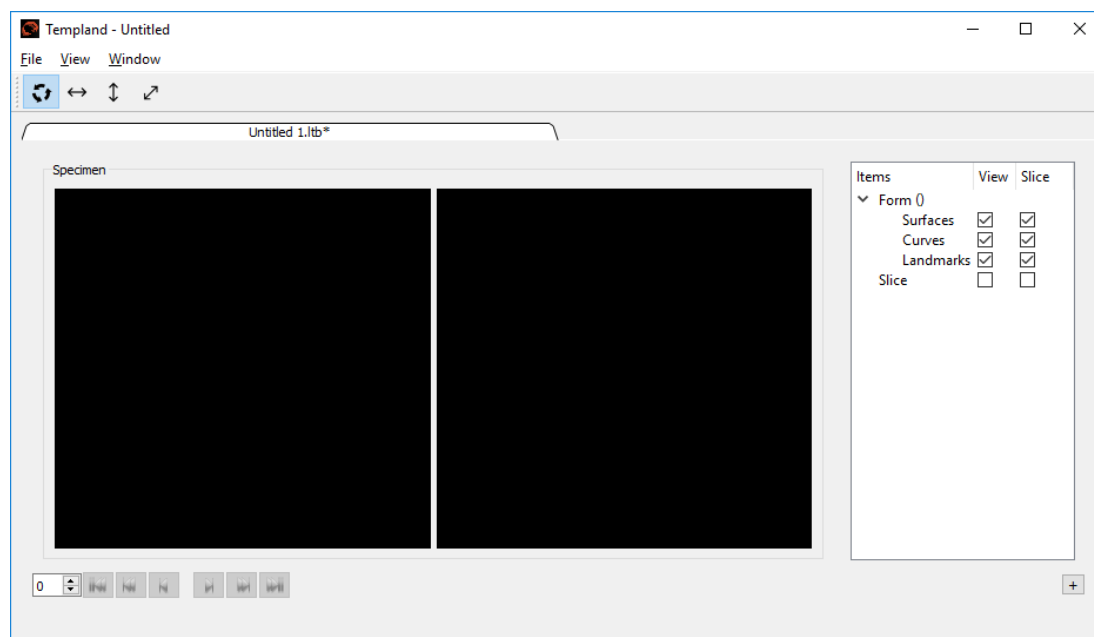

The Templand GUI is the window within which you can identify landmarks and semilandmarks. Its **left** side will show the **specimen in 3D**, its **right** side will show a particular **2D slice** through the specimen. It can be beneficial for many applications to minimize the right window. The Templand window has a number of menued items at the top.

### File Menu

The file menu allows you to:

- create a **“New”** tableau,
- **“Open File”** an existing tableau file,
- **“Save”** and **“Save as”** tableau files,
- and to **“Refresh Output Ports”**. It is important to note this item because it updates the data for the export, i.e. the LM and sLM coordinates that you have gathered in Templand and that you now want to submit to the ET Core for further analysis.

### View Menu

The view menu contains two items

- **“Template”**: toggles between the view for one specimen (as you have seen above) and a **view for two specimens**. The one **on top must be the template** specimen which is required to transfer landmarks to the target (which you must load in the bottom windows).
- **“Background”**: lets you change the background colour.

### Window Menu

The window menu allows to arrange views conveniently, as known from other Windows applications.

### The Item Menu

To the right of the two windows showing the specimen is the Item Menu. When data have been imported, this contains information about the specimen in terms of **the surface, the curves, the landmarks and the semilandmarks**. All together is called a **“form”** and can be saved as a form file (\*.frm).




The interface also allows the user to control whether or not the slice through the object is visible in either the left hand window (the view window) or the right hand window (the slice window). Surfaces, curves and landmarks can also be turned off and on by ticking the appropriate check boxes.

## 8. Loading a Surface or a Curve

The first step in building a form is to load a surface. To do this:

1. Right click on “**Surfaces**” in the Item Menu and select “Import Surface”
2. Select the surface that you wish to import from your directory (**the only surface format supported is .obj**) and choose “Open”. Then provide an “Item ID”. This name will be used in the Item Menu and **has to be consistently used for all surfaces in your sample** (for our example data, we use “Surface\_Dentin” as Item ID). The 3D surface will appear in the left hand window (View) and a section through this surface will appear in the right (Slice).
3. In the same way you can load a curve (**the only curve format supported is .obj**). The curve has to be a polyline and may be created in Templand or in another software. Right click on “**Curves**”, select “Import Curve”, and select the curve object to load. Again provide an “Item ID”. Each name will be used in the Item Menu and **has to be consistently used in your sample for a particular type of curve** (for our example data, we use “Curve\_Talon” and “Curve\_Trigon” as Item IDs). The curve(s) will show up in the left and right hand windows as a red line.

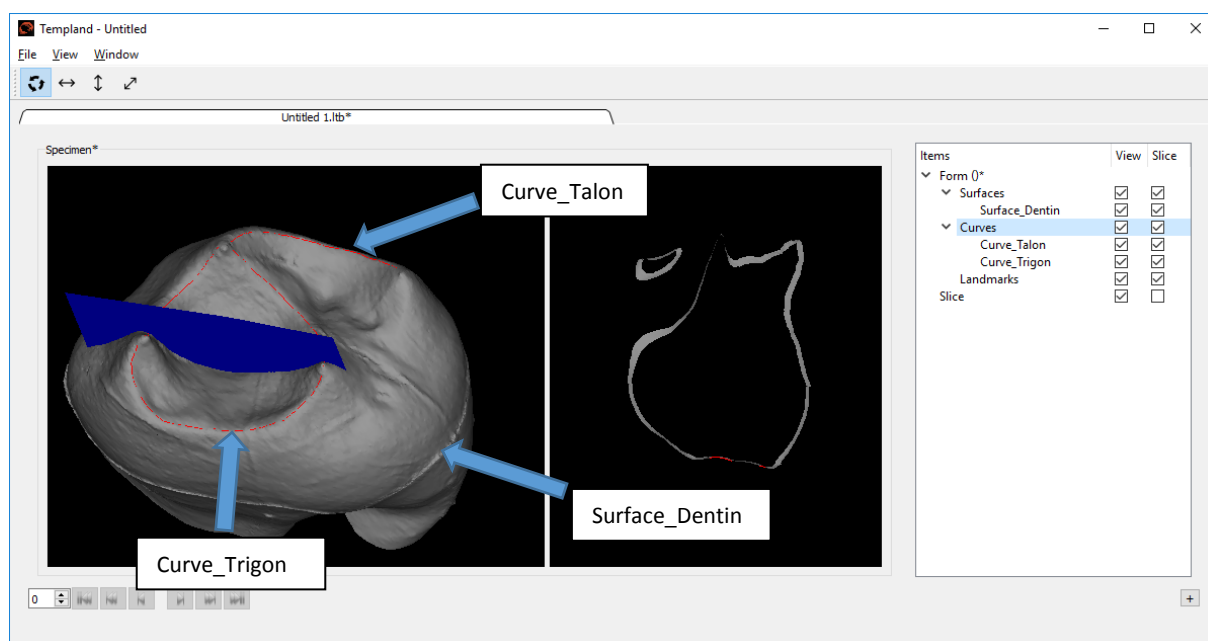

Activating the “View” for “Slice” in the Item Menu will show the actual slice position for the right hand window in the 3D rendering of the left hand window as a blue plane. In both windows the view can be modified following the instructions below:

### Templand User Controls

#### Mouse

- |                              |                     |
|------------------------------|---------------------|
| Left mouse button drag       | ⇒ Tumble            |
| Left mouse button drag + Alt | ⇒ Rotate            |
| Left mouse double click      | ⇒ Create a Landmark |

Right mouse button drag	⇒ Zoom
Right mouse button	⇒ Selects a Landmark
Middle mouse button drag	⇒ Translate
Middle mouse button scroll	⇒ Clipping
<b>Keys</b>	
Arrow Keys	⇒ 90 degree rotations
r R	⇒ 120 degree rotation
[ ]	⇒ 90 degree twists
+ -	⇒ Jump Zoom
Page Up/Down	⇒ Jump Clipping
M	⇒ Move selected Landmark
C	⇒ Center View
Spacebar	⇒ Focus at View
<b>Pressing and holding Ctrl</b> along with keys/mouse navigation buttons will result in <b>fine control</b> of the movement.	

**Note:** the direction and keyboard keys only work when surface or slice is highlighted in the items menu; it will not work if landmarks are highlighted. Use the jump push function when zooming in on a position the camera view cuts through the surface you wish to magnify. In the Templand GUI top right underneath the menu labels there are a **series of arrows** which can be selected. These constrain the way in which the slice moves under mouse control.



## 9. Placing a landmark

To **place a landmark**, position the cursor over the relevant position in either the view (left) or the slice (right) window and **double click**. Check the point from different viewing angles if it is correct, it might look different from different views. A **red cross** will appear at the location of the landmark and in the **Item Menu**, the landmark will be **named** starting with LM00001, LM00002, ..., LM99999. If you want to give it a proper name (which can make sense, e.g. for classic landmarks), then right click on the landmark's current name and choose "Rename" and type a new name. Left of the name appears **[F]** which indicates that the position of this landmark is "**fixed**" which would be the correct status for a classic landmark located by the user. It was not "**warped**" **[W]** from another specimen, and it was not "**projected**" **[P]** or "**slid**" **[S]**. You will see these attributes later in the context of semilandmarks.

If there is more than one landmark, you can select each landmark in the Item Menu. A currently selected landmark or the most recently placed landmark is displayed in red, the remaining landmarks in green. The colours and the symbols that are used to display the landmarks can be changed by right clicking on the relevant landmarks tab within the Item Menu and selecting "Colour" or "Symbol" (cross, circle, or dot). The size cannot be changed in the current version.

If you want to **delete a landmark**, right click on it in the Item Menu and select "Delete". You should do this only for the last landmark you located because you will otherwise **mess up the order and naming of your landmarks!** In case you are not satisfied with the position of a landmark that isn't last, right click on it in the Item Menu and select instead "**Move**". Drag your mouse to the new position and double click. The landmark will be moved there.

## 10. Adding semilandmarks

If you wish to work with semilandmarks in addition to classic landmarks, **for one specimen** in your sample – the template – **you have to locate those semilandmarks**. They can be on curves and on surfaces and will be arranged in groups. In order to add semilandmarks, right click on the “Landmarks” heading within the Item Menu and select “**Add Semilandmark**”. You must now provide an “**Item ID**” for the first group of semilandmarks. This name will be used in the Item Menu and **has to be consistently used in your sample for a particular type of semilandmarks** (for our example data, e.g. we use “sLM\_Curve\_Talon”, “sLM\_Curve\_Trigon”, and “sLM\_Surface\_Trigon” as Item IDs). Most importantly, by choosing from a list, **embed** (affiliate) your semilandmark group **with a particular item** which has to be an object already loaded in your form. In our example data

“sLM_Curve_Talon”	are embedded in	“Curve_Talon”
“sLM_Curve_Trigon”	are embedded in	“Curve_Trigon”
“sLM_Surface_Trigon”	are embedded in	“Surface_Dentin”.

Click on a semilandmark group in the Item Menu to highlight it. Now you can start placing semilandmarks along a curve or on a surface the same way you located landmarks. The automated numbering will start with SM00001, SM00002, ..., SM99999. You can rename semilandmarks though this doesn’t make much sense since semilandmarks are treated as a group. Once all semilandmarks for one semilandmark group have been defined, you can create another group. In our example we will use two curves and one surface with sLM.

## 11. Save and export a form file

At any stage during landmarking you can **save the form file**, which will contain the coordinates of each individual landmark and semilandmark, the pointer to the objects used (surfaces, curves), and the embedding information. Right click on “Form” (top line in the Item Menu) and select “Save” or “Save as”.


The **coordinates** of landmarks and semilandmarks can be **exported** at any time as simple text file which is easy to read with any software, and of course also ET Core. Right click on “Landmarks” (Item Menu) and select “**Export Landmarks for this Specimen**”. The text file will have no labels, just X-, Y-, and Z-coordinates per landmark and semilandmark, but the order of LM and sLM will be exactly that one of the Item Menu.

## 12. Creating a tableau file

A tableau is a sequence of form files. It can be used to export data of the sample or even to show the landmarking process. It can be created at any point in time. However, it is advisable to **save a tableau at the beginning of a process**, e.g. before sliding all targets against the template (see chapter “5. The principal workflow” above). To create a tableau file select from the top menu bar “File” and “Save as” to save the new tableau file (.ltb). In the case multiple forms have been appended to the tableau file, their coordinates can be exported in one process. Right click on “Landmarks” (Item Menu) and select “**Export Landmarks for all Specimens**” to **export X- Y-, and Z-coordinates of all landmarks and semilandmarks for all specimens** included in the tableau file.

**ET Core can also read tableau files directly**, thus an export is only necessary if you intend to use the coordinate data in another software. To import a tableau file, create an “Import” node in ET Core, double click it and **select “Tableau” as “Datafile format”**. This will import all coordinate data of all specimens but not any surface file. As

usual (see “Manuals for the EVAN Toolbox No. 2”), the surface of one particular specimen has to be loaded separately. ET Core will also accept .stl beside .obj for visualization.

The tableau file can be replayed as it is being built (or reloaded from a saved tableau file) by using the controls at the bottom left of the Templand window. These **controls**  are **similar** to the controls you would find them on an **MP3-player** and have similar action. The tableau allows also to record the process of landmarking:

1. Load a surface and curves
2. Prior to adding any landmarks select “File” and “Save as” to save the new tableau file (.ltb)
3. Add one landmark as described above
4. Right click on “Form” in the Item Menu and select “Append”
5. Add a second landmark
6. Right click on “Form” in the Item Menu and select “Append”
7. Following each subsequent landmark added select “Append” before adding the next landmark
8. After the landmark sequence is complete “Save” the tableau file in its current stage

**Note:** each individual landmark will be saved as a .frm file, and there will also be a tableau file save in the .ltb format. It is the tableau (.ltb) file that needs to be opened in order to replay the landmark sequence.

### **13. Working with the example data**

As outlined above, the surface files for the dentin crowns of three human upper second molars have to be created in other software. Our sample includes “152\_uLM2\_Dentin.obj” which will be the *template*, “697\_uLM2\_Dentin.obj” as *target no. 1*, and “S4\_uLM2\_Dentin.obj” as *target no. 2*.

In the following instructions, **input files** will be printed in **blue**, **output files** in **green**.

As next step we describe how to create a curve in ET Templand.

#### **Create a curve**


Open Templand

Right click “Surfaces” from the Item Menu

Select “Import Surface” and select

“152\_uLM2\_Dentin.obj” as file > for “Item ID” specify “Surface\_Dentin”

Locate landmarks for the curve as described above under “Landmarking”. The curve in ET will be a polyline which is connecting points with straight lines (it is not interpolating a spline). The more points you define, the better the representation will be (see image to the right). In our example, place approximately 70 LM on the trigon ridge, start at the Metacone, move clockwise, and always stay on the ridge.

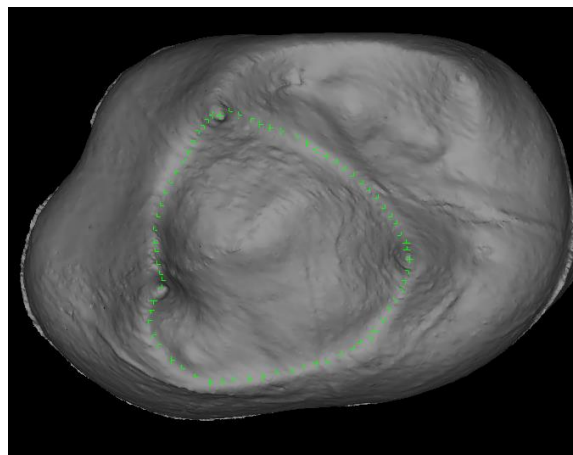

Right click “Landmarks” > “Save as Curve” > name it “152\_Curve\_Trigon.obj”

Once the curve is saved, you do not need these 70 landmarks anymore.

Delete all the landmarks that were defining the polyline by right clicking “Landmarks” > “Delete all Landmarks”

Create in the same way another curve for the talon which you call “152\_Curve\_Talon.obj”.

Create the curves for the other specimens in the sample.

### **Locate classic landmarks on the template**

Now that we obtained polylines for curves, the actual process of landmarking starts.

We have to define classic landmarks prior to the placing of semilandmarks on curves and surfaces. This step is obligatory, since LMs are required to make the TPS warping of semilandmarks possible.

Open Templand

Right click "Surfaces" from the Item Menu

Select "Import Surface" and select "[152\\_uLM2\\_Dentin.obj](#)" as file > as "Item ID" specify "Surface\_Dentin"

Click on "Landmarks" in the Item Menu

Locate the 10 classic landmarks on the object as outlined in "3. The example data" and rename them accordingly.

Right click "Form" > "Save as" > name it "[152dentin\\_trueLM.frm](#)"

Study the structure of the .frm file in Wordpad, you might need to interact with it directly in case of problems.

### **Import curves for the template**

We can now import the predefined curves to prepare our form file for the adding of semilandmarks.

Right click "Curves" > "Import Curve" > select "[152\\_Curve\\_Talon.obj](#)" > as "Item ID" specify "Curve\_Talon"  
The curve will appear on the surface as a red line. You can change colour by right-clicking "Curve\_Talon" > "Colour".

Right click "Curves" > "Import Curve" > select "[152\\_Curve\\_Trigon.obj](#)" > as "Item ID" specify "Curve\_Trigon"

Right click "Form" > "Save as" > name it "[152dentin\\_trueLM&2curves.frm](#)"

### **Add semilandmarks to the template**

As described in "10. Adding semilandmarks" we have to locate them once on the template specimen. Our goal here is to include semilandmarks on two curves and on the trigon surface.

Right click "Landmarks" > "Add semilandmarks" > as "Item ID" specify "sLM\_Curve\_Talon" > Embed with item > "Curve\_Talon"

Select "sLM\_Curve\_Talon" in the Item Menu

Place semilandmarks along the polyline of the talon as outlined in "3. The example data" (green)

*If you are not exactly on the polyline when you place them, the semilandmarks will be automatically projected to the nearest point of the polyline (you will recognize the [P] in front of the semilandmark group). If you are not satisfied with the location of a semilandmark, select it from the Item Menu, hold "Ctrl+M" keys, and move your mouse over a different location. Click the left mouse button once to place it there.*

Right click "Landmarks" > "Add semilandmarks" > as "Item ID" specify "sLM\_Curve\_Trigon" > Embed with item > "Curve\_Trigon"

Select "sLM\_Curve\_Trigon" in the Item Menu

Place semilandmarks along the polyline of the trigon as outlined in "3. The example data" (blue)


Right click "Landmarks" > "Add semilandmarks" > as "Item ID" specify "sLM\_Surface\_Trigon" > Embed with item > "Surface\_Dentin"

Select "sLM\_Surface\_Trigon" in the Item Menu

Place semilandmarks in the basin of the trigon as outlined in "3. The example data" (red)

*When double-clicking the left mouse button, the surface semilandmarks will be automatically projected to the nearest point of the surface. In this case they keep the attribute [P] because they will be projected and slid later for all target specimens. If you are not satisfied with the location of a placed semilandmark, select it from the Item Menu, hold "Ctrl+M" keys, and move your mouse over that semilandmark. It will follow now your mouse pointer. Click the left mouse button once to place it in a different location.*

Right click "Form" > "Save as" > name it "[152dentin\\_trueLM&2curves&surface\\_sLM.frm](#)"

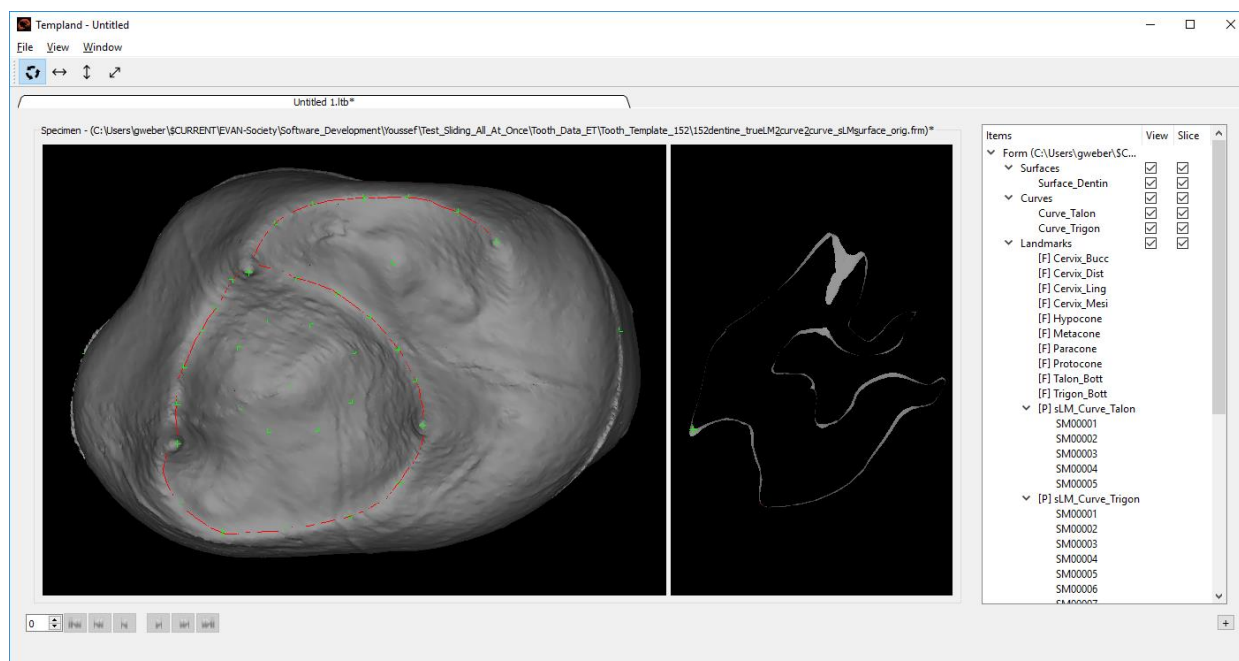

### Preparing the target specimens

The target specimens “697” and “S4” are other two upper second molars that will be treated in a similar way than the template “152”. The difference is, however, that **only the classic landmarks have to be defined** and the predefined **curves have to be imported** into the form files. The semilandmarks will then be warped automatically from the template specimen.

#### Open Templand

Right click “Surfaces” from the Item Menu

Select “Import Surface” and select “697\_uLM2\_Dentin.obj” as file > as “Item ID” specify “Surface\_Dentin”

Click on “Landmarks” in the Item Menu

Locate the 10 classic landmarks on the object as outlined in “3. The example data” and rename them accordingly.

Right click “Form” > “Save as” > name it “697dentin\_trueLM.frm”

Right click “Curves” > “Import Curve” > select “697\_Curve\_Talon.obj” > as “Item ID” specify “Curve\_Talon”

Right click “Curves” > “Import Curve” > select “697\_Curve\_Trigon.obj” > as “Item ID” specify “Curve\_Trigon”

Right click “Form” > “Save as” > name it “697dentin\_trueLM&2curves.frm”

Right click “Form” > “Reset”

Right click “Surfaces” from the Item Menu

Select “Import Surface” and select “S4\_uLM2\_Dentin.obj” as file > as “Item ID” specify “Surface\_Dentin”

Click on “Landmarks” in the Item Menu

Locate the 10 classic landmarks on the object as outlined in “3. The example data” and rename them accordingly.

Right click “Form” > “Save as” > name it “S4dentin\_trueLM.frm”

Right click “Curves” > “Import Curve” > select “S4\_Curve\_Talon.obj” > as “Item ID” specify “Curve\_Talon”

Right click “Curves” > “Import Curve” > select “S4\_Curve\_Trigon.obj” > as “Item ID” specify “Curve\_Trigon”

Right click “Form” > “Save as” > name it “S4dentin\_trueLM&2curves.frm”

### **Create sliding semilandmarks of the target specimens**

All the form files for the three specimens in our sample are now ready for the last and crucial step, the warping of the semilandmarks from the template to the targets and subsequent iterative slide-project step.

Open Templand

Click "File" > "New"

Right click "Form" > "Load" > select file "[152dentin\\_trueLM&2curves&surface\\_sLM.frm](#)"

Click "File" > "Save as" > type "[uLM2.Itb](#)" to save a new tableau file

Right click "Form" > "Append" to append the template as the first specimens to the tableau

Click "View" > "Template" to change to the 4-window view

In the upper row, right click "Form" > "Load" > select file

"[152dentin\\_trueLM&2curves&surface\\_sLM.frm](#)"

In the lower row, right click "Form" > "Reset"

In the lower row, right click "Form" > "Load" > select file

"[697dentin\\_trueLM&2curves.frm](#)"


In the upper row, right click "Landmarks" > "Warp all unmatched landmarks"

*A dialogue window will appear that allows you to specify the parameters for the slide-project step. Activating the checkbox for "Slide Semi-Landmarks" enables the process. The options are to "Slide-Project Once", or to specify the number of iterations under "Slide-Project # iterations", or to define a certain difference from one iteration to the next using "Slide-Project until error".*

Enable sliding by checking the box on top


Choose one of the options, e.g. Slide-Project # iterations = 7

Click "OK"

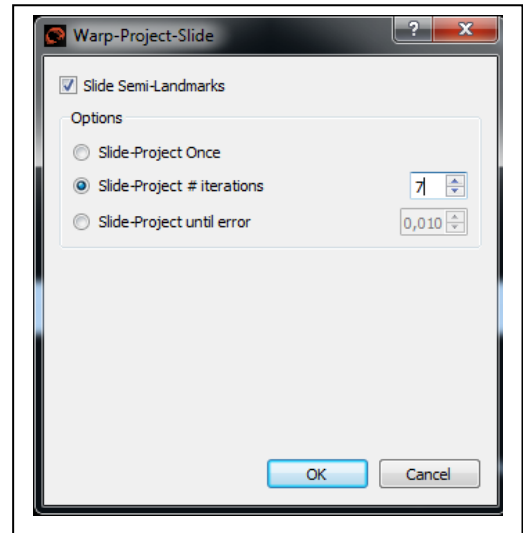

*The unmatched semilandmarks, i.e. "sLM\_Curve\_Talon", "sLM\_Curve\_Trigon", and "sLM\_Surface\_Trigon" are now warped to the target specimen and then slid and projected 7 times. Depending on the speed of the processor and the number of landmarks this may take several seconds to minutes. The result appears in the lower window. The target specimen now has all semilandmarks in the Item Menu and they are shown on the surface of the specimen in the View window.*



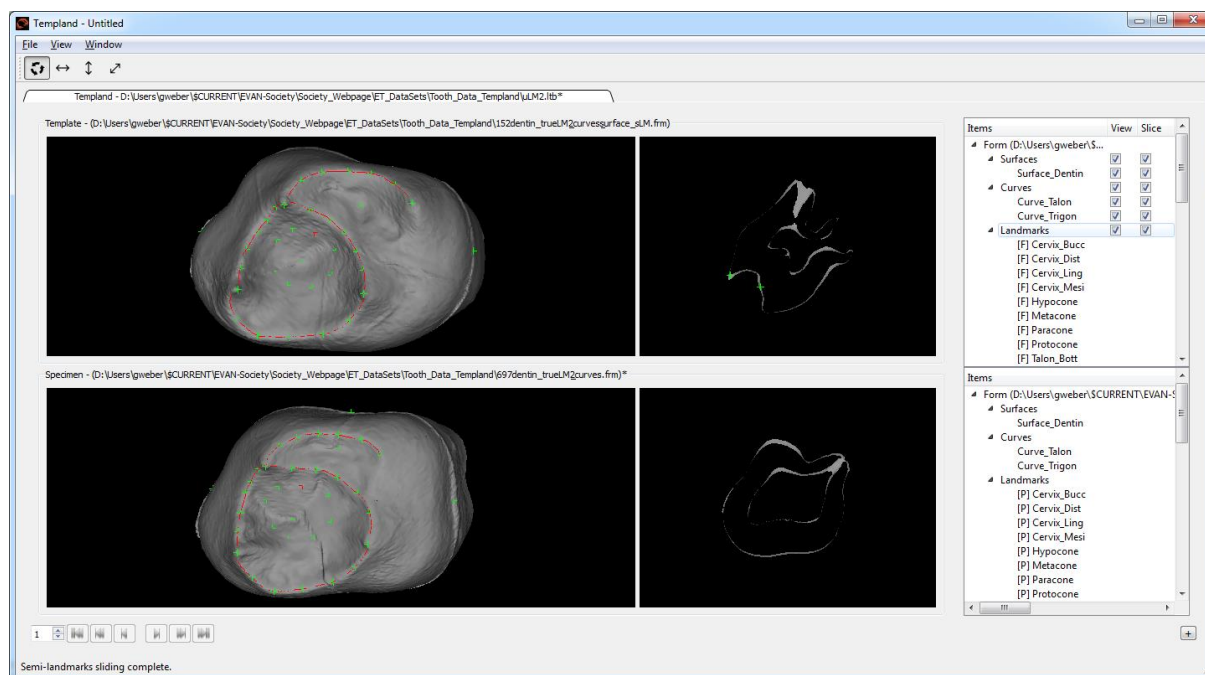

In the lower row, right click "Form" > "Append" to append the target no. 1 as the second specimen in the tableau

Click "File" > "Save" to save the tableau in its current state

In the lower row, right click "Form" > "Reset"

In the lower row, right click "Form" > "Load" > select file "[S4dentin\\_trueLM&2curves.frm](#)"

In the upper row, right click "Landmarks" > "Warp all unmatched landmarks"

Enable sliding by checking the box on top

Choose one of the options, e.g. Slide-Project # iterations = 7

Click "OK"

In the lower row, right click "Form" > "Append" to append the target no. 2 as the third specimen in the tableau

Click "File" > "Save" to save the tableau in its current state

### **Export the Landmark configurations for analyses**

There are several ways of getting data from ET Templand out to ET Core or to other programmes.

1. A simple way to store the created coordinate data is **in a text-file**:



In the lower row, right click on "Landmarks" > Export Landmarks for all Specimens" > name it "[uLM2.txt](#)"


```
'D:\path\152dentin_trueLM&2curves&surface_sLM.frm
4.27058 11.0306 3.92942
9.56317 10.3978 10.6095
12.4628 6.34926 10.8249
8.80493 6.28765 3.49794
7.92563 6.13324 12.3458
... .
```

This file can be read with any editor and be easily imported by other spreadsheet or statistic applications. The order of landmarks corresponds exactly with the order in the Item Menu of Templand.

2. Coordinates can be taken over by ET Core and be exported there **in different formats, e.g. Excel or Morphologika**



Click "File" > "Refresh Output Ports"  
*(this can take a while, depending on the number of LM and performance of the computer)*  
 Switch back to ET Core window  
 Place an "Export" node on the workspace  
 Connect "All Specimens" of the "Templand" node with "Specimens" of the "Export" node  
 Run the network   
 Double click "Export" node > Specimens "Morphologika" format > "Export Specimens" > provide file name "uLM2\_morphologika.txt" > OK  
 Pause the network 

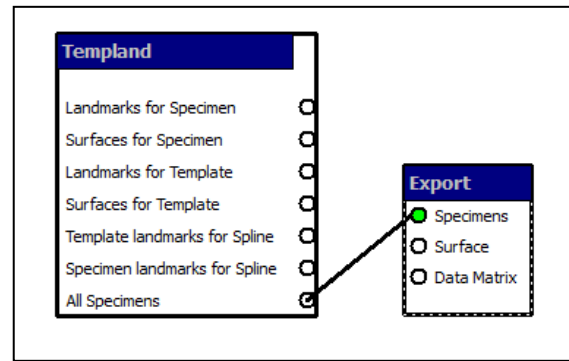

When you open the textfile, you can read and edit the data file in Morphologika format:

```
[individuals]
3
[landmarks]
38
[dimensions]
3
[names]
1
2
3
[rawpoints]
' 1
4.27058141475 -3.92942137468 11.030601127
9.56316928506 -10.6095383808 10.3978187865
12.4628461435 -10.8248545249 6.34925518826
8.80493200822 -3.49794476756 6.28765330736
7.92562694985 -12.3458311982 6.13323854725
...
```

3. Data from the **tableau file can be imported directly in ET Core** at any point in time.


If you not yet familiar with ET Core, please consider the instructions in "Manuals for the EVAN Toolbox No. 2" for the use of "Import" and "GPA".

Use a predefined "GPA.vpn"

Double click the "Import" node and select "Tableau" as "Datafile format" > select file "uLM2.ltb"

Run the network 

Open the "3DViewer" to get a visualization of the three upper second molars after a Generalized Procrustes Analysis (see below). You can also import one of the surfaces that were used for this example. From here, you can proceed with all kinds of shape and form analyses in ET Core.

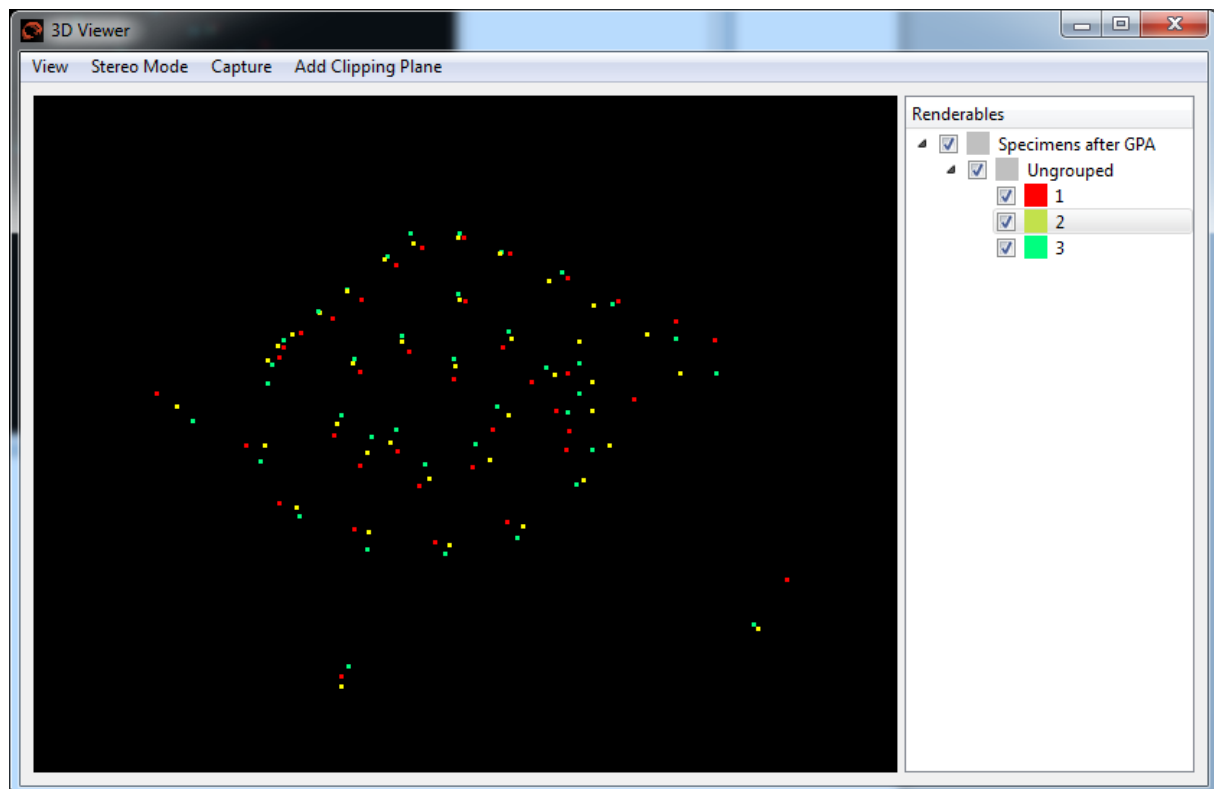

Production of the EVAN Toolbox and preparation of this manual were supported by Marie Curie Research Training Network EVAN (MRTN-CT-2005-019564) and the EVAN-Society.

Copyright by the EVAN-Society.